


# Ambu® IV Trainer

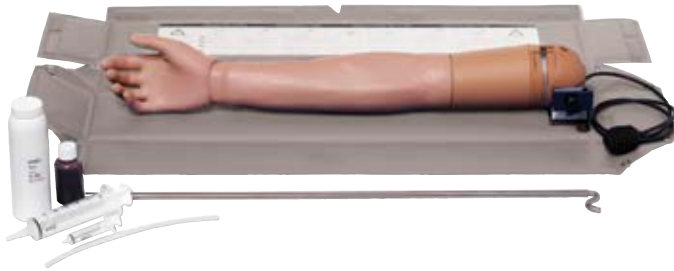

## Key features

- Replaceable skin and veins
- Movable wrist and rotating arm
- Indicator panel
- Multi purpose carrying case
- Fits all AmbuMan type manikins suitable for the connection of arms
- Valve system

## Ambu® IV Trainer

The Ambu IV Trainer is for training the insertion of cannulas and catheters, the infusion of fluids and the injection of medication, blood sampling and pulse measurement.

### Replaceable skin and veins

Replaceable skin and veins are made of natural latex, which features hole closure after penetration, extending the lifetime of the arm.

### Movable wrist

The movable wrist and rotating arm ensure an absolutely realistic training product, especially when a catheter is placed on the back of the hand.

### Indicator panel

The indicator panel, which represents each blood vessel, gives immediate feedback, e.g. if an infusion has been correctly accomplished.

### Valve system

The valve system reduces the formation of air pockets in blood vessels.

## Specifications


Weight	3 kg
Length	68 cm
Reservoir capacity	260 ml
Carry case dimensions	75 x 15 x 11 cm

## Order information

Item no.	Description
255 001 000	Ambu IV Trainer, complete with infusion stand, red dye concentrate, talcum powder and 5 ml and 50 ml syringes. Supplied in carry case / washable underlay
<b>Accessories</b>	
255 000 704	Red dye concentrate

## Materials

Product	Material
Skin	Natural latex
Arm and hand	PVC, hard
Veins and artery	Natural latex
Indicator	Silicone
Storage bag	Nylon

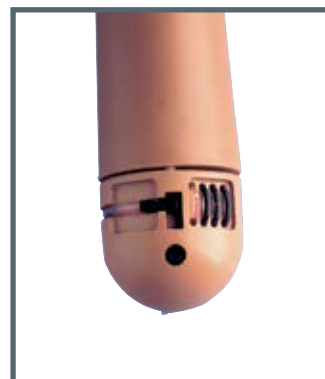

Indicator panel!



Accessories for Ambu IV Trainer



Packed in a robust bag

### Ambu A/S

Baltorpbakken 13  
DK - 2750 Ballerup  
Denmark  
Tel. +45 7225 2000  
Fax +45 7225 2053  
www.ambu.com